
# PULSONIC

powered by **BRAUN**

S 26.523.3
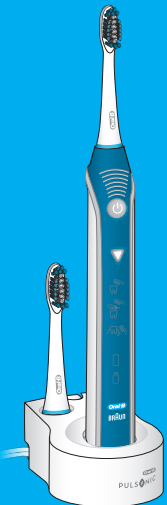
|               |    |       |  |
|---------------|----|-------|--|
| Deutsch       | 4  | DE AT | 00 800 27 28 64 63                                       |
| English       | 7  | CH    | 08 44 - 88 40 10   |
| Français      | 10 | GB    | 0 800 783 7010   |
| Español       | 14 | IE    | 1 800 509 448  |
| Português     | 17 | FR    | 0 800 944 803<br>(appel gratuit depuis un<br>poste fixe) |
| Italiano      | 20 |       |  |
| Nederlands    | 23 | BE    | 0 800 14 592   |
| Dansk / Norsk | 26 | LU    | 800 21172  |
| Svenska       | 29 | ES    | 901 11 61 84   |
| Suomi         | 32 | PT    | 808 20 00 33   |
| Ελληνικά      | 35 | IT    | 800 124 600  |
|               |    | NL    | 0800-4 45 53 88  |
|               |    | DK    | 70 15 00 13  |
|               |    | NO    | 22 63 00 93  |
|               |    | SE    | 020 - 21 33 21   |
|               |    | FI    | 020 377 877  |
|               |    | GR    | 800 801 3457   |

**Internet:**

[www.oralb.com](http://www.oralb.com)  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)

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CH: 0844 - 88 40 10

# English

**Oral-B® Pulsonic** has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, some basic safety precautions should always be followed.

## IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age three. Children from 3 to 14 years, persons with reduced physical, sensory or mental capabilities, or persons with lack of experience or knowledge of the product can use the toothbrush, but should have supervision and/or instruction concerning use of the appliance in a safe way, and understand the hazards involved.
- Cleaning and maintenance should not be performed by children.
- Children should not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

## WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.

- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

## Description

- a Oral-B Pulsonic brush head
- b On/off button
- c Brushing mode button
- d Handle
- e Brushing mode indicator
- f Charge indicator light
- g Low charge indicator light
- h Charging unit (basic charger and charger cover with release button)

## Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit into an electrical outlet and place the toothbrush handle on the charging unit (see picture ①, page 3). The charge indicator light (f) flashes while the handle is being charged. Once it is fully charged, the charge indicator light shines continuously for 5 seconds; then it turns off. A full charge takes at least 12 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- When the rechargeable battery is low, the red low charge indicator light (g) starts flashing after turning off the toothbrush. You will also hear a sound reminding you to charge your toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. There is no risk of overcharging the battery. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

For electric specifications, please refer to the bottom of the charging unit.

## Brushing

### Brushing modes

Your toothbrush offers different brushing modes for varying oral care needs:

- «Daily Clean» – Exceptional mouth cleaning.
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas.
- «Massage» – Gentle massage of gums.

- When pressing the on/off button (b), your toothbrush automatically starts in the «Daily Clean» mode.
- To switch to the other modes, successively press the brushing mode button (c).
- Your toothbrush memorizes the brushing mode even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the brushing mode resets to the «Daily Clean» mode.
- To turn your toothbrush off, press the on/off button (see picture ②, page 3).

The brushing mode is displayed in the Brushing mode indicator display (e) while brushing.

### Brushing technique

- Slide the brush head on the handle. Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on your toothbrush.
- Place the brush head on the outside surface of your upper teeth. The toothbrush bristles should be placed against the teeth at a slight angle towards the gumline.
- Turn your toothbrush on and start brushing in a slightly circular motion. Apply light pressure during brushing. After a few seconds, guide the bristles to the next section. Brush inside, outside and chewing surfaces of your teeth with the same motion throughout the mouth for the recommended 2 minutes (see picture ③, page 3).
- To avoid splashing, turn your toothbrush off before removing from the mouth by pressing the on/off button.

You may experience a slight tickling or tingling sensation when you first use your Oral-B Pulsonic. This will subside as your teeth and gums will become accustomed to your Oral-B Pulsonic cleaning experience. If you experience this initial tingling, you may want to brush in the «Sensitive» mode first.

In the first few days of using your toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

### Brush heads

#### Oral-B Pulsonic brush head

The Oral-B Pulsonic brush head (a) is recommended for everyday deep cleaning. Its bristles provide a thorough cleaning between your teeth, tooth surfaces and along the gumline. The blue INDICATOR bristles help you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approxi-

mately within 3 months, signalling it is time to replace your brush head.

### Professional Timer

The «Professional Timer» helps you brushing all four quadrants of your mouth equally. It signals with short beep sounds at 30-second intervals that it is time to move to the next quadrant of your mouth. A single beep indicates completion of the first quadrant, two beeps indicate the second quadrant, three beeps indicate the third quadrant and four beeps sound when you reach the recommended 2-minute mark (see picture ④, page 3).

Your «Professional Timer» memorizes the elapsed brushing time, even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds, the timer resets to the initial starting point.

### Cleaning recommendations

- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Wipe them dry before you put the handle on the charging unit.
- From time to time you can also clean the charging unit with a damp cloth. For thorough cleaning you may remove the charger cover (h) by pushing the release button on top of the charging unit. The charger cover is dishwasher-safe.
- The charger unit should never be placed in water (see picture ⑤, page 3).

Subject to change without notice.

### Environmental Notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Center or at approved recycling or disposing locations in accordance with local regulations.

### Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or

replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

UK:  
FREEPOST RSHJ-KRSA-UHGE  
Braun Product Service  
Olympic Shaver Centre Ltd  
PO BOX 9799  
Calverton  
Nottingham  
NG14 6WN

ROI:  
Braun Product Service  
Olympic Shaver Centre Ltd  
Dept AA1758  
PO BOX 6151  
FREEPOST  
Dublin 2

## 30-Days Money Back Guarantee

Try Oral-B Pulsonic starting from the day of purchase. If you are not satisfied with this product, please return your handle, Refills and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

Great Britain / Ireland:  
«Oral-B 30-Day Challenge»  
Dept. PX415  
PO Box 10  
Diss  
Norfolk  
IP98 3SP

For further details, please contact our customer service.

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