




Register your product and get support at
www.philips.com/welcome

SCH540

PHILIPS
AVENT

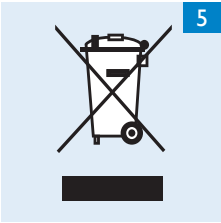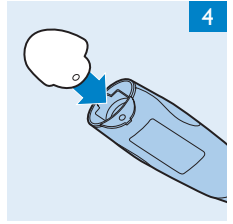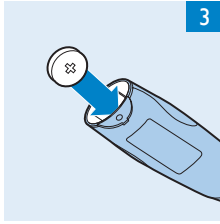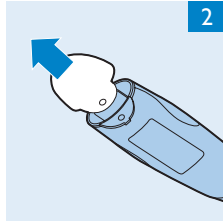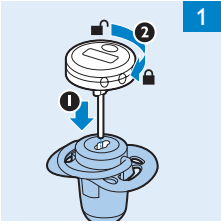ENGLISH	6
DANSK	13
DEUTSCH	20
ΕΛΛΗΝΙΚΑ	28
ESPAÑOL	37
SUOMI	45
FRANÇAIS	52
ITALIANO	60
NEDERLANDS	68
NORSK	76
PORTUGUÊS	83
SVENSKA	91
TÜRKÇE	98

Introduction

Congratulations on your purchase and welcome to Philips AVENT! To fully benefit from the support that Philips AVENT offers, register your product at www.philips.com/welcome.

Philips AVENT is dedicated to producing caring, reliable products that give parents the reassurance they need.

The Philips AVENT Digital Thermometer Set accurately and conveniently measures your baby's body temperature orally, axially and rectally. The regular soother and the soother thermometer are suitable for use with children aged 0-12 months. The digital flexible-tip thermometer is suitable for children and adults of all ages.

The Digital Thermometer Set is designed for your baby's safety, comfort, well-being and, therefore, your own peace of mind. The unique combination of a soother thermometer and a digital, flexible-tip thermometer offers you the choice of adapting the temperature measurement method to suit your baby.

The storage box contains three items:

- A Philips AVENT soother thermometer to reliably measure your baby's temperature. It is very convenient to use and comfortable for your baby.
- A regular Philips AVENT Freeflow soother to let your baby get used to the shape. You can use this soother as a regular soother.
- A Philips AVENT digital flexible-tip thermometer to measure your baby's body temperature axially or rectally if you want additional reassurance.

Make sure you follow the safety guidelines when using these items so that you do not expose your child to any risk.

Use the soother thermometer to obtain a temperature indication. Use the flexible-tip thermometer to determine the actual body temperature. If in any doubt, consult your doctor.

Important

Read this user manual carefully before you use the thermometers and save it for future reference.

Danger

- Never attach ribbons or cords to a soother; your child may be strangled by them.
- Keep the protection caps for the teats away from children to avoid suffocation.
- Inspect the soothers carefully before each use, especially when the child has teeth. Pull the soothers in all directions. Throw away at the first signs of damage or weakness.

Warning

- Do not immerse the measurement unit of the soother thermometer in water or any other liquid. Never sterilise or boil the measurement unit.
- Never leave your child when you are measuring its temperature with the soother thermometer or the flexible-tip thermometer.
- The thermometers are not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the thermometers by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the thermometers.
- Never use the soother thermometer measurement unit separately to avoid injury.
- Battery compartment covers must be replaced properly if removed. If loose, they may choke your child.

Caution

- We recommend replacing the regular soother after 4 weeks, for safety and hygiene reasons.
- Always sterilise or boil the teat, shield and protection cap of the soother and of the soother thermometer after use (see chapter 'Cleaning').

- Do not leave the soother, soother thermometer and the flexible-tip thermometer in direct sunlight or near a source of heat.
- Do not leave the soothers in disinfectant ('sterilising solution') for longer than recommended, as this may weaken the teat.
- Do not drop the thermometers or expose them to shocks.
- Do not bend or disassemble the measurement unit of the soother thermometer or the tip of the flexible-tip thermometer
- Do not use the soother thermometer as a regular soother as this reduces its useful life.
- Do not use the thermometer if the display shows the battery symbol. This means that the battery is low and the measurements will be inaccurate.
- Only use these thermometers to measure human body temperatures.
- Never dip the teats in sweet substances or medication, your child may get tooth decay.
- In the event that either soother becomes lodged in the mouth, **DO NOT PANIC**. The soothers cannot be swallowed and are designed to cope with such an event. Remove from the mouth with care, as gently as possible.
- When not in use, store the soother, the soother thermometer and the flexible-tip thermometer in the storage box supplied.
- Do not use this Philips AVENT soother thermometer in combination with parts of the old Philips soother thermometer (SCH530).

Electromagnetic fields (EMF)

These Philips AVENT thermometers comply with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the thermometers are safe to use based on scientific evidence available today.

Before first use

Clean the soother, soother thermometer and the flexible-tip thermometer before first use. See chapter 'Cleaning'.

Using the appliance

Ways to use the thermometers

The most convenient way to measure your baby's temperature is to use the soother thermometer according to this user manual. This method gives a reliable indication. If you need additional reassurance (for example if the first measured temperature is above 38°C), you can then use the digital flexible-tip thermometer to take an axillary (underarm) or a rectal (anal) measurement.

Note: Oral measurements are usually 0.5 to 1°C lower than the actual body temperature.

Childrens' temperature

The temperature of a child is normally between 36.5 and 37.5°C. However, you have to keep in mind that a child's temperature varies depending on age, activity, time of day, etc. Also remember that some illnesses do not have high temperature as a symptom. Temperature measurement should not be a substitute for consulting your physician. When reporting your child's temperature to a professional, state the measurement method.

Soother thermometer

The most convenient way to measure your baby's temperature is to use the soother thermometer. This is also the most comfortable method for the baby.

Do not use the soother thermometer for 30 minutes after the baby has had hot or cold drinks, as this causes an inaccurate measurement.

Do not use the soother thermometer as a regular soother.

- 1 Check if the measurement unit of the soother thermometer is firmly attached to the soother (i.e. turned completely to the closed lock position).

If the measurement unit is not firmly attached to the soother, the measurement may be inaccurate and the measurement unit may become a suffocation hazard.

- 2 Place the soother in the baby's mouth for approx. 3 minutes to let the thermometer warm up.

If you do not do this, the measured temperature may be inaccurate.

- 3 Press the on/off button  to switch on the thermometer.

Note: You can either press the on/off button while the soother is still in your baby's mouth, or you can remove the soother, press the on/off button and place the soother back into your baby's mouth.

- 4 Leave the soother in the baby's mouth until the thermometer beeps.

See section 'Audible signals' at the end of this chapter for the meaning of the beeps.

- ▶ The display now shows the measured temperature and the °C indication lights up continuously.

Note: If for some reason the temperature cannot be measured, the display shows an 'L' or an 'H'. This means that the temperature is outside the temperature range. Switch off the thermometer, check if it is correctly positioned and switch it on again to restart the measurement.

- 5 Remove the soother thermometer from the baby's mouth.

Do not leave the child alone with the soother thermometer.

- 6 Press the on/off button  to switch off the thermometer.

Note: Switch off the thermometer after use to save battery power. If you do not switch off the thermometer, it switches off automatically after 10 minutes.

Digital flexible-tip thermometer

Do not insert the flexible-tip thermometer into the ear or nose.

The soother thermometer is just as accurate as the digital flexible-tip thermometer and is more convenient to use. However, because rectal measurement is inherently more accurate than oral measurement, the digital flexible-tip thermometer gives you an opportunity to reassure yourself. For a rectal measurement of a baby's temperature, proceed as follows:

- 1 Remove the protection cap from the digital flexible-tip thermometer.

Make sure that the metal probe of the flexible tip is still firmly attached to the thermometer.

Tip: If available, use special plastic hygienic covers to cover the tip of the thermometer (these are available from most chemists).

- 2 Lubricate the tip of the thermometer (or hygienic cover) with a water-soluble jelly.

Water-soluble jelly is available from most chemists.

Note: Do not use petroleum jelly.

- 3 Press the on/off button  to switch on the thermometer.

- 4 Lie the baby on its stomach or its side.

- 5 With one hand, separate the baby's buttocks. With the other hand, gently slide the tip of the thermometer no more than 12mm (1/2 inch) into the rectum. If you detect any resistance, stop immediately.

Hold the baby firmly so that it cannot turn while the thermometer is inserted.

Never push the flexible tip too far into the anus.

6 Keep the thermometer in place until it beeps.

See section 'Audible signals' at the end of this chapter for the meaning of the beeps.

Note: If for some reason the temperature cannot be measured, the display shows an 'L' or an 'H'. This means that the temperature is outside the temperature range. Switch off the thermometer, check if it is correctly positioned and switch it on again to restart the measurement.

7 Remove the thermometer and read the display.**8** Press the on/off button  to switch off the thermometer.

Note: Switch off the thermometer after use to save battery power. If you do not switch off the thermometer, it switches off automatically after 10 minutes.

Measuring the temperature of older children

For older children and adults, the digital flexible-tip thermometer can also be used to measure temperature orally (in the mouth) or axially (in the armpit).

Do not use the flexible-tip thermometer for oral temperature measurements on babies.

If the thermometer has been used rectally or in the armpit without hygienic covers, do not use it for oral measurement anymore for hygienic reasons, not even after cleaning.

Audible signals

Thermometer	Signal	Meaning
Soother thermometer	4 beeps	Measured temperature is normal (below 37.5°C)
	3 sets of 4 short beeps	Possible fever (measured temperature above 37.5°C)
Flexible-tip thermometer	10 beeps in 10 seconds	Measured temperature is normal (below 37.5°C)
	3 sets of short beeps in 10 seconds	Possible fever (measured temperature above 37.5°C)

Features**Last temperature recall**

Both the soother thermometer and the digital flexible-tip thermometer have a memory function. If you want to compare temperatures, you can recall the previously measured temperature. This can only be done once before a new measurement.

1 Press and hold the on/off button  for about 3 seconds.

▶ The display shows the last measured temperature and an 'M' for memory.

Note: When you release the on/off button, the last measured temperature is erased.

Battery-low indication

Both the soother thermometer and the digital, flexible-tip thermometer have a long-life battery which should last for about 2 to 3 years. If the display shows a battery symbol, the battery is low. Replace the battery of the flexible-tip thermometer (see chapter 'Replacing the battery').

Note: The battery of the soother thermometer cannot be replaced. If the battery is empty, discard the soother thermometer.

Automatic shut-off

Both the soother thermometer and the digital flexible-tip thermometer switch off automatically if the thermometer has not been used for 10 minutes.

Error function

Both the soother thermometer and the digital flexible-tip thermometer have an error function. In case of malfunction, the display shows an 'E'. Proceed as follows if this happens:

- 1 Press the on/off button  to switch off the thermometer.
- 2 Press the on/off button  again to switch on the thermometer again.

Cleaning

Do not use abrasive cleaning agents, anti-bacterial cleaning agents or chemical solvents. Excessive concentration of cleaning agents may eventually cause plastic components to crack. Do not use cracked parts.

Clean the regular soother, the soother thermometer and the flexible-tip thermometer after each use to ensure optimal hygiene.

Be careful not to damage the parts when you clean them.

Cleaning the soother thermometer

Before first use, wash and rinse the soother thoroughly.

Clean the soother before each use.

Before cleaning, always remove the measurement unit from the teat and shield.

- 1 Turn the measurement unit 90° anticlockwise with respect to the shield.
- 2 Pull the measurement unit straight out of the teat and shield.

Note: Be careful not to bend the tip of the measurement probe.

- 3 Clean the thermometer measurement unit with a soft damp cloth.

Do not immerse the measurement unit in water.

Do not sterilise or boil the measurement unit.

Do not clean the measurement unit in the dishwasher, in the microwave or in a microwave steriliser.

- 4 Sterilise or boil the soother and protection cap (not the measurement unit!).

The soother and cover are suitable for all forms of sterilisation.

Sterilise using an AVENT Steam Steriliser or boil in water for 5 minutes.

Do not use abrasive cleaning agents or anti-bacterial cleaners. Excessive concentration of detergents may eventually cause plastic components to crack. Should this occur, replace immediately.

If you use a sterilising solution, follow the instructions on the package. Do not leave the soother in the sterilising solution for longer than recommended, as this may weaken the teat.

Wash your hands thoroughly and ensure surfaces are clean before contact with sterilised components.

- 5 After sterilisation, let the soother cool down for 5 minutes and squeeze out any trapped water from the teat. This is to ensure hygiene.

To reassemble the soother thermometer:

- 1 Align the blocks of the measurement unit with the matching slots in the back of the shield.

- 2 Insert the measurement unit into the shield. Make sure that the tip of the probe goes into the pocket inside the teat. (Fig. 1)
- 3 Turn the measurement unit 90° clockwise to lock it into position.

Cleaning the flexible-tip thermometer

Do not sterilise or boil the thermometer.

Do not clean the flexible-tip thermometer in the dishwasher, in the microwave or in a microwave steriliser.

- 1 Clean the flexible-tip thermometer with a soft damp cloth.

Cleaning the regular soother

Before first use, wash and rinse the soother thoroughly.
Clean the soother before each use.

- 1 **Sterilise or boil the soother.**

The soother and cover are suitable for all forms of sterilisation.

Sterilise using an AVENT Steam Steriliser or boil in water for 5 minutes.

Do not use abrasive cleaning agents or anti-bacterial cleaners. Excessive concentration of detergents may eventually cause plastic components to crack. Should this occur, replace immediately.

If you use a sterilising solution, follow the instructions on the package. Do not leave the soother in the sterilising solution for longer than recommended, as this may weaken the teat.

Wash your hands thoroughly and ensure surfaces are clean before contact with sterilised components.

- 2 **After sterilisation, let the soother cool down for 5 minutes and squeeze out any trapped water from the teat. This is to ensure hygiene.**

Cleaning the storage box

Do not sterilise the storage box.

- 1 Wipe the storage box with a damp cloth.

Storage

You can use the storage box to store the soother, the soother thermometer and the flexible-tip thermometer or to take them with you when you travel.

- 1 Put the protection caps on the soother, the soother thermometer and the flexible-tip thermometer after use to keep them clean.
- 2 Put the soother, the soother thermometer and the flexible-tip thermometer in the storage box.
- 3 Put the storage box in a safe, dry place away from direct sunlight and heat sources.

Replacing the battery

Swallowing a button battery is extremely dangerous.

Batteries may explode if exposed to high temperatures or fire.

Never attempt to recharge the battery.

The thermometers have a long-life battery that should last about 2 to 3 years. If the display shows a battery symbol, the battery is low.

Replace the battery of the flexible-tip thermometer as described below.

Note: The battery of the soother thermometer cannot be replaced. If the battery is empty, discard the soother thermometer.

- 1 Remove the battery compartment cover (Fig. 2).
- 2 Use an insulated probe to remove the battery (e.g. a toothpick).
- 3 Insert a new battery (LR41) with the + pole facing upwards (Fig. 3).
- 4 Reattach the battery compartment cover (Fig. 4).
- 5 Dispose of the empty battery in the proper way (see chapter 'Environment').

Environment

- Non-rechargeable batteries contain substances that may pollute the environment. Do not throw away non-rechargeable batteries with the normal household waste, but dispose of them at an official collection point for batteries. Always remove the non-rechargeable batteries before you discard and hand in the thermometers at an official collection point.
- Do not throw away the thermometers with the normal household waste at the end of their life, but hand them in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 5).

Guarantee & service

If you need service or information or if you have a problem, please visit the Philips AVENT website at www.philips.com/AVENT or contact the Philips Consumer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Consumer Care Centre in your country, go to your local Philips dealer.

Technical specifications

	Digital flexible-tip thermometer	Soother thermometer	AVENT soother
Battery	1.55VDC (LR41)	Battery not replaceable	-
Battery life	> 100 operating hours	> 100 operating hours	-
Measurement range	32°C to 44°C	32°C to 44°C	-
Accuracy	± 0.1°C for 35°C to 38°C range	± 0.1°C for 35°C to 38°C range	-
Teat size	-	Suitable for newborn babies	Suitable for newborn babies
Compliance	Conforms to MDD93/42/EEC	Conforms to MDD93/42/EEC	Conforms to EN1400

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Philips Consumer Lifestyle B.V