
AQUACARE WATERFLOSSER

powered by **BRAUN**

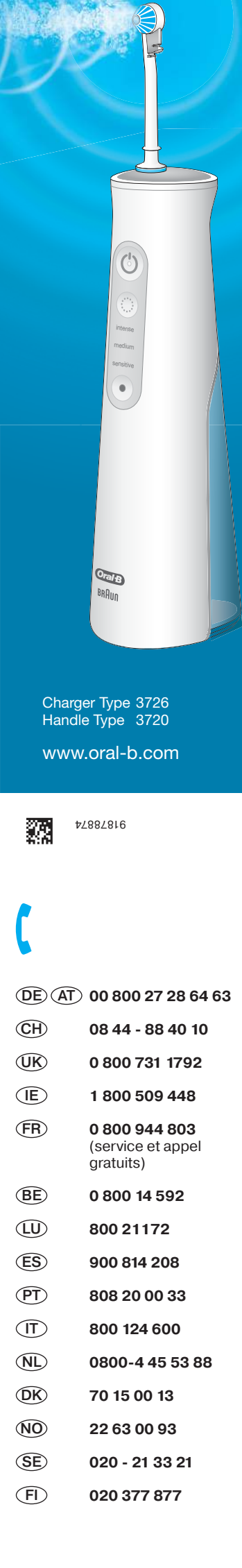Charger Type 3726

Handle Type 3720

www.oral-b.com

91878874

(DE) (AT) 00 800 27 28 64 63

(CH) 08 44 - 88 40 10

(UK) 0 800 731 1792

(IE) 1 800 509 448

(FR) 0 800 944 803
(service et appel gratuits)

(BE) 0 800 14 592

(LU) 800 21172

(ES) 900 814 208

(PT) 808 20 00 33

(IT) 800 124 600

(NL) 0800-4 45 53 88

(DK) 70 15 00 13

(NO) 22 63 00 93

(SE) 020 - 21 33 21

(FI) 020 377 877

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

www.az-oralb.it

Braun GmbH

Frankfurter Straße 145











61476 Kronberg/Germany

Charger Type 3726

Handle Type 3720

91878874/X-18

DE/EN/FR/ES/PT/IT/NL/DK/NO/SE/FI

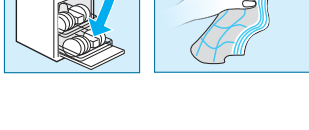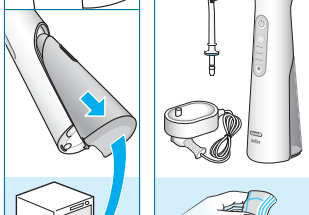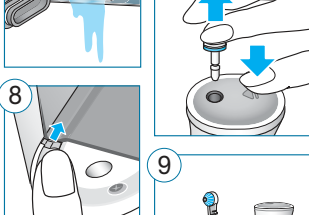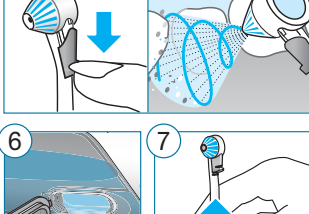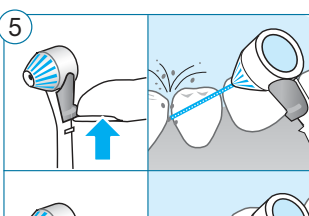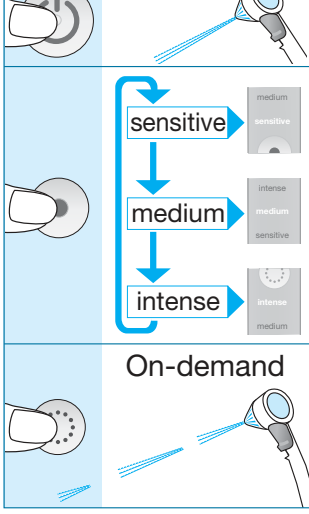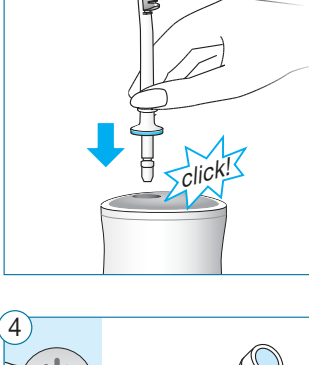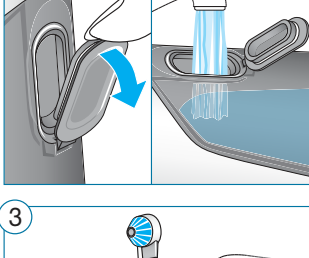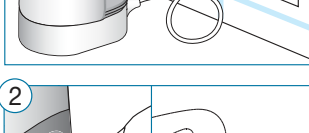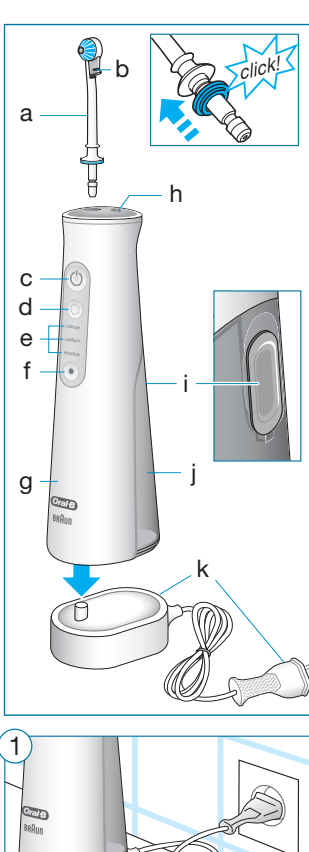Deutsch

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English

Welcome to Oral-B!

Please read the usage instructions carefully before first use and save this manual for future reference.

WARNING

- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable. Do not disassemble the product except when disposing of the battery. When taking out the battery from the handle for disposal, use caution not to short the positive (+) and negative (–) terminals.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This appliance is a personal care device and is not intended for use on multiple persons.
- Small parts may come off, keep out of reach of children.

To avoid irrigator nozzle breakage which may generate a small parts choking hazard or damage teeth:

- If the product is dropped, the irrigator nozzle should be replaced before the next use even if no damage is visible.
- Clean the appliance properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the appliance.

IMPORTANT

- Periodically check the entire product/cord for damage. If the product/cord is damaged, take it to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- The oral irrigator is not intended for use by children under age 8. Children aged from 8 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the product, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be performed by children.
- Oral irrigators shall be used only with liquids defined in the manufacturer's instructions.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in the appliance. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Only use the charger provided with your appliance. Do not charge any other device with this charger, as this may cause damage.

Description

- a Nozzle
- b Switch
- c On/off button (continuous jet)
- d On-demand button

- e Intensity level
- f Intensity selection button
- g Handle
- h Release button
- i Water reservoir cap
- j Water reservoir (150 ml)
- k Charger

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Noise level: 68 dB (A)

Operating and charging

Plug the charger (k) into an electrical outlet and place the handle (g) on the charger (picture 1).

While charging the on/off button (c) light is fading in and out.

Once it is fully charged the light turns off. A full charge typically takes 14 hours.

If the battery is running low the on/off button is flashing in red.

First Use

Attach the nozzle (a) to the handle until it locks into place indicated by a click sound (picture 3). Fill the water reservoir (j) with some fresh and clean water by lifting the reservoir cap (i) on the back of the handle (picture 2).

Hold the handle over a washbasin, press the on/off button (c) and allow the entire water volume to run through the nozzle. When the device has not been used for more than one week repeat this step to ensure flushing.

Jet functions:

You can choose between different jet stream types by sliding the jet switch (b) down or up (before doing so, turn off the water stream with the on/off button (c) on the handle (g)):

Rotating stream type for general cleaning and gum massage (jet switch in lower position): The rotation of the water stream with its micro-bubbles are designed for easy penetration and gentle gum massage (picture 5).

Focused stream type (jet switch in upper position): The focused water stream with its micro-bubbles are designed for targeted cleaning (picture 5).

Press the on/off button (c) for a continuous jet stream.

Press the on-demand button (d) to trigger a single jet stream. Upon releasing the on-demand button the single jet stream stops.

Press the intensity selection button (f) successively to adjust the intensity level (sensitive – regular/medium – intense) (depending on model). You may pre-select the intensity level before you turn on the handle. The intensity level selected last is memorized (picture 4).

Note: In the beginning, we recommend to use the sensitive level for your teeth and gums to get accustomed to the usage.

Using the irrigator

Fill the water reservoir (j) with lukewarm water by lifting the reservoir cap (i) on the back of the handle (g).

Lean over the sink and direct the nozzle in your mouth. Close your lips slightly and turn on the handle by pressing the on/off button (c).

Guide the jet along the gumline and pause briefly between teeth.

Alternatively, you may also point the nozzle tip at a desired area and use the on-demand button (d) for a single jet stream.

Keep your mouth slightly open to allow the water to run into the sink.

Refill water reservoir as many times as needed to complete your cleaning routine. It may take multiple refills to clean your entire mouth.

After use

Make sure the water reservoir is completely empty. Press the on/off button (c) and leave it turned on until no more water is running through the nozzle, to prevent any deposit build-up or bacteria from collecting in the standing water. Turn off the handle by pressing the on/off button again once the water is empty. A longer operation of the handle without water may damage the motor.

To drain any water out of the handle hold it over the washbasin and tilt it from side to side (picture 6).

Release the nozzle (a) by pressing the release button (h) and pulling the nozzle at the same time (picture 7).

Lift up the cap (i) and let the water reservoir air-dry. Always wipe the handle and nozzle dry.

Note: Make sure you always empty the water reservoir before storing or travelling with the irrigator.

In general we recommend to replace the nozzle every 3 months.

Cleaning recommendations


Clean and dry handle after each use and leave the water reservoir cap (i) open.

The water reservoir (j) is dishwasher-safe and can be detached via the snap-fit on the bottom of the handle. Move the lever and slide the reservoir down (picture 8). Regularly clean the charger and handle with a damp cloth. Store the appliance in a cool and dry place.

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

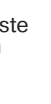Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads/nozzles as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

UK only:

This guarantee in no way affects your rights under statutory law.

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