SPCMG-06SB
SPCMG-06SR

DE	Gebrauchsanweisung	Massage-Pistole	Seite 6
FR	Mode d'emploi	Pistolet de massage	Page 12
IT	Manuale d'uso	Pistola per massaggi	Pagina 18
GB	Instructions for use	Massage Gun	Page 24

GARANTIE

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
WARRANTY


This appliance has a two year warranty period. During this period the manufacturer is responsible for any failures that are the direct result of construction failure. When these failures occur the appliance will be repaired or replaced if necessary. The warranty will not be valid when the damage to the appliance is caused by wrong use, not following the instructions or repairs executed by a third party. The guarantee is issued with the original till receipt. All parts, which are subject to wear, are excluded from the warranty.

If your device breaks down within the 2-year warranty period, you can return the device together with your receipt to the shop where you purchased it.

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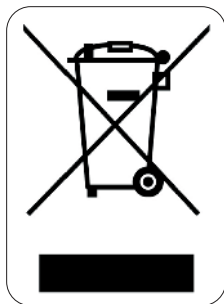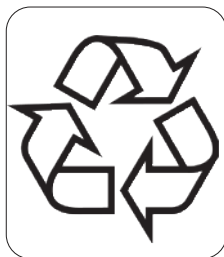RECYCLING INFORMATION

This symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it must be brought to the applicable collection point for recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. For more detailed information about recycling of this product, please contact your local city office, your household waste disposal service or the shop where you purchased the product.

The packaging is recyclable. Please treat package ecologically.

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Dear customer

We wish you much pleasure with this SPC quality product and thank you for your confidence.

The SPC brand has stood for consistency, quality and development power for almost 20 years. Each device has been carefully manufactured and tested.

With the proverbial Swiss thoroughness, accuracy and honesty, the company from St. Gallen has convinced numerous Swiss customers with innovative products since its foundation in 2001.

With over 1000 products, SPC offers a comprehensive range in the various product groups from the electrical household sector and consumer electronics. The right product for every occasion!

Get to know the variety of SPC products and visit our online shop at www.spc.ch.

Your SPC Electronics AG

Service hotline:

Switzerland:	071 227 99 99
Monday to Friday:	08.00 - 12.00h / 13.30 - 16.30h

Spare parts can be ordered at: www.spc.ch

SPC Electronics AG
Mövenstrasse 12
CH-9015 St. Gallen
www.spc.ch

SAFETY INSTRUCTIONS

When using electrical appliances, basic safety precautions should always be taken, including the following:

- Read all instructions carefully. Keep this manual for future reference.
- Make sure all packaging materials and promotional stickers have been removed before using the appliance for the first time. Make sure children cannot play with the packaging materials.
- This appliance is intended to be used in household and similar applications such as:
 - staff kitchen areas in shops, offices and other working environments;
 - farm houses;
 - by clients in hotels, motels and other residential type environments;
 - bed and breakfast type environments.
- Children should be supervised to ensure that they do not play with the appliance.
- This appliance can be used by children aged from 16 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are older than 16 and supervised.
- Keep the appliance and its cord out of reach of children younger than 16 years.
- To reduce risks of breakdown, electric shock, injury, fire, death and equipment damage, please read all instructions and warnings.
- Keep away from liquids and heat sources. Keep ventilation ports free from dust and debris.
- Unplug the charger from device battery after charging.
- Use only the applicators that provide the best desired effect.
- Use this product in accordance with the instructions.
- Check the device carefully before each use.
- For adult only. Do not use if injured. Muscle injuries can be used for recovery. Please consult your doctor before using this device.



- Please do not use the massage gun, or any vibration device without first obtaining approval from your doctor if any of the following apply:
Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device, but we advise you to consult a doctor first.
- Use only on dry and clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use this device only on the body's soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
- Do not use in areas where bruising, pain or uncomfortable part.
- Do not place any objects into the ventilation holes.
- Do not immerse in water or allow water to enter the ventilation holes.
- Do not drop or misuse the device.
- Only recharge with supplied charger.
- Do not tamper or alter the massage gun.

Do not touch this area when starting the machine without massage head assembled.

PARTS

1. Interchangeable applicators
2. Silica gel covered handle
3. Control panel
4. Gear indicator light
5. Button
6. Power indicator light
7. Round head
8. Listic head
9. U-shape head
10. Conical head
11. Flat head
12. Bullet head

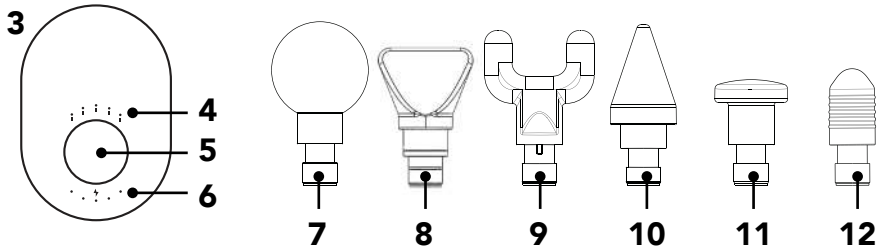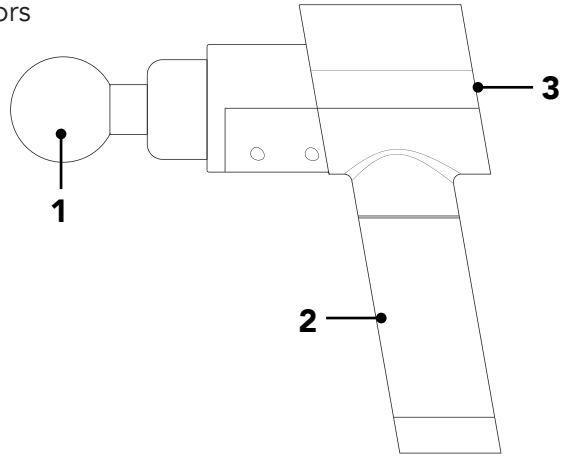USE

1. Install one of the massage heads. Rotate the massage head gently while pulling/loading, do not pull/load directly (except u-shaped massage head)
 - Round head: suitable for large muscle groups
 - Listic head: suitable for the lower back
 - U-shape head: suitable for neck & spine
 - Conical head: suitable for precise massage of pain points with high irritation
 - Flat head: suitable for any part of the body
 - Bullet head: suitable for joints, deep tissue and trigger points
2. Long press the power button once to turn device on to level one, each one more press of the button will increase the level position by 1. When it's at level 5, press the button to enter standby mode (level 0), from working mode level 5 to level 0. To restart, wait for

3 seconds and press the button to enter level 1, and repeat the previous steps.

3. In any gear, press long to enter standby mode.
4. In standby mode, long press can turn it off directly, and it will be off automatically if there is no operation for 15 seconds.

CHARGING

1. Fully charge battery for up to 6 hours before first use either attached or detached from the device.
2. To charge, connect DC end of the supplied 24V adapter to the battery's charging port and plug AC side cable into wall outlet.
3. The LED light on the device flashes regularly, indicating the battery level, indicating that it is charging.
 - When the first indicator is on, the remaining capacity of the battery is 20% to 35%;
 - When the second indicator light is on, the remaining capacity of the battery is 40% to 60%;
 - When the third indicator light is on, the remaining capacity of the battery is 65- 80%;
 - When the forth indicator light is on, the remaining capacity of the battery is 85-100%.
4. Battery may be recharged at any time and at any battery level. Do not switch on to use this device while charging.
5. The battery needs to be recharged when it runs out to less than 5% (the first LED flashes).
6. The average run-time is 4+hours depending on speed level and pressure applied during use.

CLEANING AND MAINTENANCE

1. Use slightly damp towel cleaning equipment, wipe dry with a soft cloth, do not use detergent, strong detergent, concentration more than 75% of alcohol, disinfectant and other irritating chemicals wipe clean.
2. For storage or travel, remove battery and store in the supplied box.