







PULSONIC

SLIM

powered by BRAUN

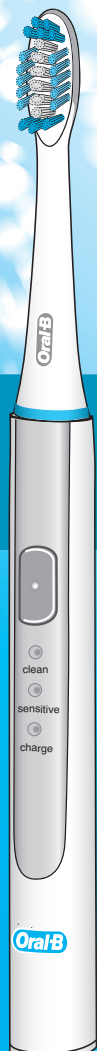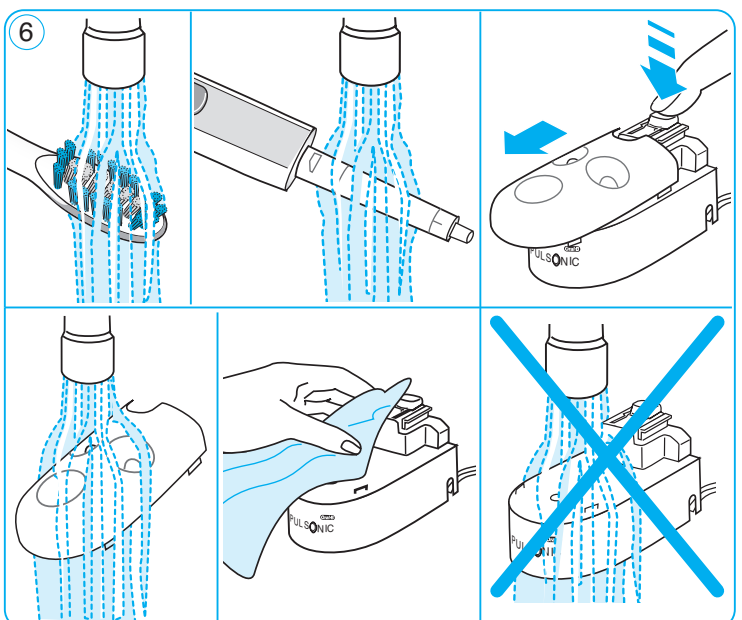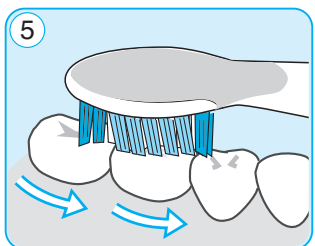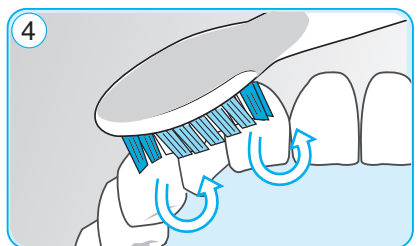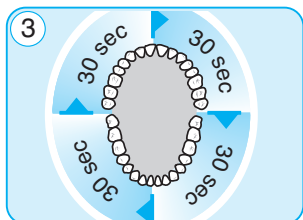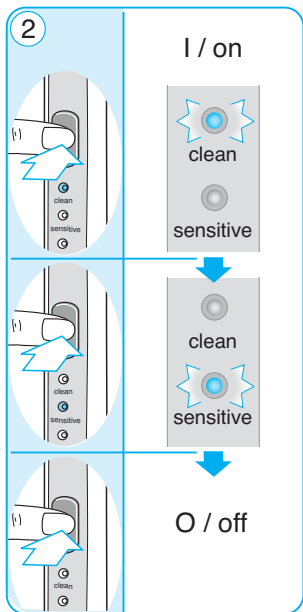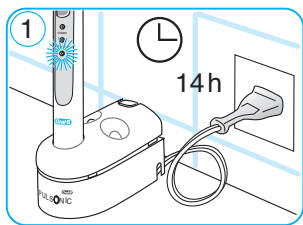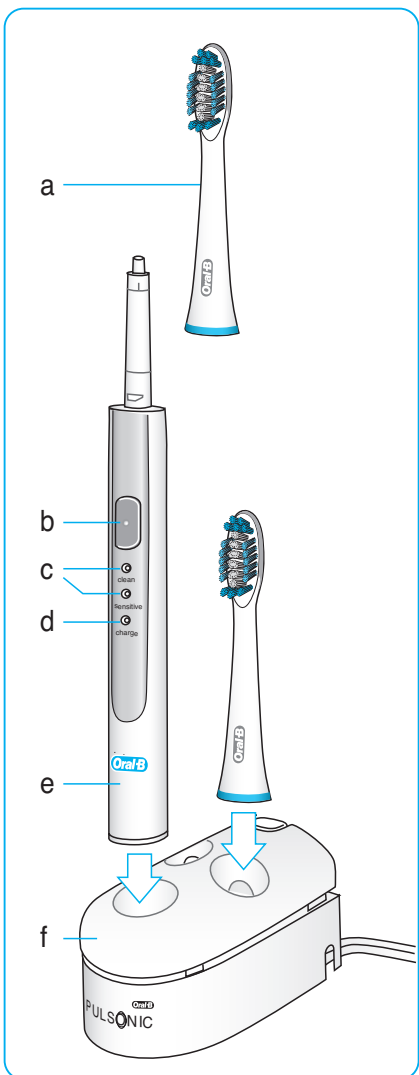Deutsch	3
English	4
Français	5
Español	7
Português	9
Italiano	10
Nederlands	11
Dansk / Norsk	13
Svenska	14
Suomi	15
Ελληνικά	17

Internet:
www.oralb.com
www.braun.com
www.service.braun.com

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CH: 0844 - 88 40 10

English

Oral-B Pulsonic Slim has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, basic safety precautions should always be followed.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral B Braun Service Center. A damaged or non-functioning unit should no longer be used.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify the product. Do not disassemble or repair. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service center. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Never drop the appliance or insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments or chargers not recommended by the manufacturer.
- Keep cord away from heated surfaces.
- Do not put metal objects such as coins and clips on the charger. This may generate heat damaging the unit and heating the metal item presenting a potential burn hazard.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

- a Oral-B Pulsonic brush head
- b On/off button
- c Mode indicators
- d Charge indicator
- e Handle
- f Charging unit with handle stand and release button

Connecting and charging

Your toothbrush is electrically safe and designed for use in the bathroom.

- Plug the charging unit (f) into an electrical outlet and place the handle (e) on the charging unit (see picture ①, page 2).
- A full charge takes at least 14 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible.

However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.

- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Brushing

Brushing modes

Oral-B Pulsonic Slim offers different brushing modes for varying oral care needs:

- «Clean» – Gentle, effective mouth cleaning.
- «Sensitive» – Gentle cleaning for sensitive areas.

When pressing the on/off button (b), your toothbrush automatically starts in the «Clean» mode. To shift from the «Clean» to the «Sensitive» mode, press the on/off button once again. To turn your toothbrush off, press the on/off button (see picture ②, page 2).

Brushing technique

- Slide the brush head on the handle. Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on your toothbrush.
- Place the brush head on the outside surface of your upper teeth. The toothbrush bristles should be placed against the teeth at a slight angle towards the gumline.
- Turn your toothbrush on and start brushing in a slightly circular motion (see picture ④, ⑤, page 2). Apply light pressure during brushing. After a few seconds, guide the bristles to the next section. Brush inside, outside and chewing surfaces of your teeth with the same motion throughout the mouth for the recommended 2 minutes.
- To avoid splashing, turn your toothbrush off before removing from the mouth by pressing the on/off button.

You may experience a slight tickling or tingling sensation when you first use your Oral-B Pulsonic Slim. This will subside as your teeth and gums will become accustomed to your Oral-B Pulsonic Slim cleaning experience. If you experience this initial tingling, you may want to brush in the «Sensitive» mode first.

In the first few days of using your toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

Brush heads

Oral-B Pulsonic brush head

The Oral-B Pulsonic brush head (a) is recommended for everyday deep cleaning. Its bristles provide a thorough cleaning between your teeth, tooth surfaces and along the gumline. The blue indicator bristles help you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling it is time to replace your brush head.

Timer

A short stuttering sound at 30-second intervals (one break) reminds you to brush all four quadrants of your mouth equally (see picture ③, page 2).

A long stuttering sound (two breaks) indicates when you reach the professionally recommended 2-minute mark.

Cleaning

Switch off the handle and remove the brush head. Rinse both brush head and handle thoroughly under running water. Wipe them dry before you put the handle on the charging unit. From time to time you can also clean the charging unit with a damp cloth. For thorough cleaning you may remove the handle stand by pushing the release button on top of the charging unit. The charging unit should never be placed in water (see picture ⑥, page 2).

Subject to change without notice.

Environmental Notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at approved recycling or disposing locations in accordance with local regulations.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

30-Days Money Back Guarantee

Try Oral-B Pulsonic Slim starting from the day of purchase.

If you are not satisfied with this product, please return your handle, Refills and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

Great Britain / Ireland: «Oral-B 30-Day Challenge»;

Braun Cashback (BC 024) PO Box 504 Leicester LE94 0AE

For further details, please contact our customer service.

This guarantee in no way affects your rights under statutory law.

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