PULSONIC

SMART SERIES

powered by BRAUN

S 32.533.5

Deutsch	4	DE AT	00 800 27 28 64 63
English	11	CH	08 44 - 88 40 10
Français	17	GB	0 800 731 1792
Español	25	IE	1 800 509 448
Português	32	FR	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	39		
Nederlands	46	BE	0 800 14 592
Dansk / Norsk	53	LU	800 21172
Svenska	60	ES	901 11 61 84
Suomi	67	PT	808 20 00 33
Ελληνικά	74	IT	800 124 600
		NL	0800 445 53 88
		DK	70 15 00 13
		NO	22 63 00 93
		SE	020 - 21 33 21
		FI	020 377 877
		GR	800 801 3457

**Internet:**

[www.oralb.com](http://www.oralb.com)  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)

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CH: 0844 - 88 40 10

## English

**Oral-B® Pulsonic SmartSeries with SmartGuide™** has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, some basic safety precautions should always be followed.

### IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age three. Children from 3 to 14 years, persons with reduced physical, sensory or mental capabilities, or persons with lack of experience or knowledge of the product can use the toothbrush, but should have supervision and/or instruction concerning use of the appliance in a safe way, and understand the hazards involved.
- Cleaning and maintenance should not be performed by children.
- Children should not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

### SmartGuide batteries

- Batteries may leak if empty or not used for a long time. Avoid skin contact whilst handling leaking batteries.
- Keep away from children. Do not swallow. If swallowed, seek medical advice at once.
- Remove fully discharged batteries immediately. Replace all batteries at the same time. Do not mix used and

new batteries, different brands or types. Do not disassemble, recharge or dispose of in fire.

### WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

### Important Information

- To avoid electromagnetic interference and/or compatibility conflicts, deactivate the radio transmission of your toothbrush before carrying or using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals (see «Deactivating the radio transmission»).
- Persons with pacemakers should always keep the toothbrush more than six inches (15 centimeters) from the pacemaker when turned on. If at any time you suspect interference is taking place, deactivate the radio transmission of your toothbrush before using it (see «Deactivating radio transmission»).

### Description


---

- a Oral-B Pulsonic brush head
- b On/off button
- c Brushing mode button
- d Handle
- e Radio transmission display
- f Charge indicator light
- g Low charge indicator light
- h Charging unit (basic charger and charger cover with release button)
- i SmartGuide display unit
- j Wall mount
- k Adhesive tape


## Quick Starter Guide

Before brushing first time, charge your toothbrush. Plug the charging unit (h) into an electrical outlet and place the handle (d) on the charging unit (see «Connecting and charging» for more details).

- Your SmartGuide display unit (i) comes with a demonstration mode activated. To leave this mode remove the battery compartment cover on the back and push the set or the h/min button, located in the battery compartment. Please note that the SmartGuide display will self-activate this demonstration mode when you do not use your toothbrush for 28 days. You can avoid this self-activation by setting the clock or synchronizing your handle.

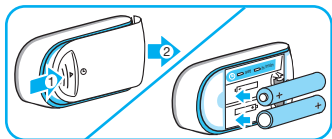- A flashing «12:00» on the SmartGuide display indicates that the time can now be adjusted. Press the «set» button, located in the battery compartment, until «12h» is displayed.

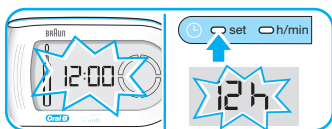- By briefly pressing the «h/min» button, you can toggle between a «12h» or «24h» daytime format. Push the «set» button to confirm your choice.

- Then the hour digit starts flashing. Press the «h/min» button until the respective hour appears and confirm with the «set» button. Proceed the same way to set both minute digits. If you prefer not using this clock feature, the flashing «12:00» will disappear automatically after 24h.

Attach the Pulsonic brush head (a) to the handle. You can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth

before switching on the appliance. Start your toothbrush by pressing the on/off button (b). To turn your brush off, press the on/off button again.

## Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit into an electrical outlet and place the toothbrush handle on the charging unit. The green charge indicator light (f) flashes while the handle is being charged. Once it is fully charged, the green charge indicator light shines continuously for 5 seconds; then it turns off. A full charge takes at least 12 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- When the rechargeable battery is low, the red low charge indicator light (g) starts flashing after turning off your toothbrush. You will also hear a sound reminding you to charge your toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. There is no risk of overcharging the battery. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

For electric specifications, please refer to the bottom of the charging unit.

## Brushing

### Brushing modes

Your toothbrush offers different brushing modes for varying oral care needs:

- |               |  |
|---------------|--|
| «Daily Clean» | – Exceptional mouth cleaning.  |
| «Sensitive»   | – Gentle, yet thorough cleaning for sensitive areas.                     |
| «Massage»     | – Gentle massage of gums.  |
| «Polish»      | – Polishing action for natural whitening for occasional or everyday use. |

«Extended Clean» – Exceptional mouth cleaning with extended brushing time. Within this mode the brushing time has been extended to 45 seconds per quadrant for a total brushing time of 3 minutes.

- When pressing the on/off button, your toothbrush automatically starts in the «Daily Clean» mode.
- To switch to the other modes, successively press the mode button (c).

- To return from any mode to the «Daily Clean» mode, press and hold the mode button.
- Your toothbrush memorises the brushing mode even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the brushing mode resets to the «Daily Clean» mode.
- To turn your toothbrush off, press the on/off button.

### Brushing technique

- Slide the brush head on the handle. Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on your toothbrush.
- Place the brush head on the outside surface of your upper teeth. The toothbrush bristles should be placed against the teeth at a slight angle towards the gumline.
- Turn your toothbrush on and start brushing in a slightly circular motion. Apply light pressure during brushing. After a few seconds, guide the bristles to the next section. Brush inside, outside and chewing surfaces of your teeth with the same motion throughout the mouth for the recommended 2 minutes.
- To avoid splashing, turn your toothbrush off before removing from the mouth by pressing the on/off button.


You may experience a slight tickling or tingling sensation when you first use your Oral-B Pulsonic SmartSeries. This will subside as your teeth and gums will become accustomed to your Oral-B Pulsonic SmartSeries cleaning experience. If you experience this initial tingling, you may want to brush in the «Sensitive» mode first.

In the first few days of using your toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

### Brush head

#### Oral-B Pulsonic brush head

The Oral-B Pulsonic brush head (a) is recommended for everyday deep cleaning. Its bristles provide a thorough cleaning between your teeth, tooth surfaces and along the gumline. The blue indicator bristles help you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling it is time to replace your brush head.

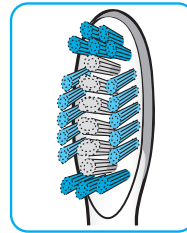
### Personalizing the timer

You can choose between the «Professional» or «2-Minute» timer and the «Count Up» or «Count Down» timer to meet your specific needs.

The «Professional Timer» helps you brushing all four quadrants of your mouth equally. It signals with short beep sound at 30-second intervals (45-second intervals within the Extended Clean mode) that it is time to move to the next quadrant of your mouth. A single beep indicates completion of the first quadrant, two beeps indicate the second quadrant, three beeps indicate the third quadrant and four beeps sound when you reach the recommended 2-minute mark. In addition the «quadrant circle» on the SmartGuide display unit helps to display your brushing progress.

The «2-Minute Timer» signals with four beeps and a «smiley» on the SmartGuide display unit

that the dentist recommended 2-minute brushing time has been reached.

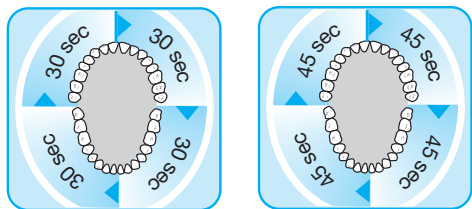
Both timers memorize the elapsed brushing time, even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds, the timer resets to the initial starting point.

The «Count up Timer» displays your actual brushing time. The «Count down Timer» indicates the remaining brushing time. Starting with «2:00» minutes, it goes down to «0:00». When brushing longer than two minutes, it automatically changes to the «Count up Timer».

The total brushing time can be extended to 3 minutes by using the «Extended Clean» mode (see «Brushing modes» for more details).

Your toothbrush comes with the «Professional Timer» and the «Count up Timer» activated. If you wish to change these settings, proceed as follows:

- Press and hold the brushing mode button (c) for 3 seconds until the SmartGuide display unit (i) shows «:30» and «timer».

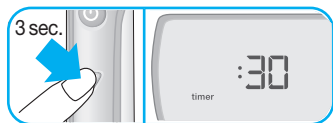- By briefly pressing the on/off button (b), switch to the «2-Minute Timer». Confirm your choice by pushing the mode button (c).

- This will lead you directly to the «Count up Timer» display («up» and «timer» will be shown). Switch to the «Count down Timer» by

pressing the on/off button and confirm it with the mode button.

- To leave the setup menu press the mode button again. The setup mode will be ended automatically if no button is pressed for 60 seconds. Please note that all 3 pilot lights are flashing while being in the setup menu.

## SmartGuide display unit

### Displayed Icon Meaning

«Daily Clean» mode

«Sensitiv» mode

«Massage» mode

«Polish» mode

«Extended Clean» mode

«Quadrant circle» – currently you are brushing in the first 30-second interval (45-second interval within Extended Clean mode). Move to the next quadrant of your mouth, when the second segment starts flashing.

Dentist recommended 2-minute brushing time (3-minute brushing time within the Extended Clean mode) has been reached.

The rechargeable battery of the handle is low. Charge the toothbrush handle.


Batteries of the SmartGuide display unit are low. Change batteries of the SmartGuide display unit.

## Wall mount

If you prefer to mount the SmartGuide display unit on a wall, mirror or other surface, use the adhesive tape (k) to fix the wall mount (j). Before doing so, make sure the chosen surface is clean and dry (remove grease and dust first).


Position the wall mount so, that the Oral-B logo is readable. Then, press the wall mount evenly and firmly to the wall. Wait 24h before sliding the SmartGuide display unit into the fixed wall mount.

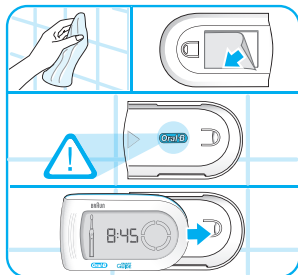
Please note: The adhesive tape will not work on dirt-repellent surfaces.

## Cleaning recommendations

- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean brush head and handle separately under running water. Wipe them dry before you put the handle back on the charging unit.

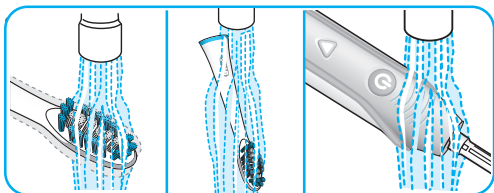- From time to time you can also clean the charging unit with a damp cloth. For thorough cleaning you may remove the charger cover (h) by pushing the release button on top of the charging unit. The charger cover is dishwasher-safe.

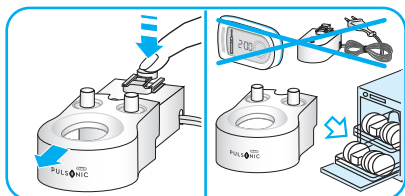
- Charging unit and SmartGuide display unit should never be placed in water.

## Synchronizing

To avoid interference of display messages when using more than one Oral-B Pulsonic SmartSeries toothbrushes with separate SmartGuide display units in the same household, you need to assign the handles to their respective SmartGuide display units.

To do so, please follow the steps below:

- Push the «h/min» button, located in the SmartGuide battery compartment, for at least 3 seconds. The SmartGuide display unit shows «L---», indicating the «Learn» mode for synchronization is activated.
- Now you are requested by a flashing on/off button in the toothbrush icon on the SmartGuide display unit to turn the handle on. When doing so, the automatic synchronization process starts. It is finished when «L-1-» is displayed.

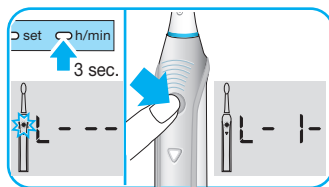
- Then turn off the handle.

If you would like to assign a second handle to the same SmartGuide display unit, repeat the last two steps described above with the other handle. You can exit the «Learn» mode by pressing the «h/min» or «set» button. Otherwise the «Learn» mode will automatically end after 30 seconds.

Please note:

- A maximum of 2 handles can be synchronized with one SmartGuide display unit. After you synchronized the second handle, «Full» will be displayed.

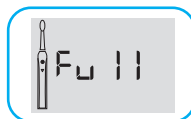- When changing batteries of the SmartGuide display unit, stored data will be erased. Therefore, please repeat the synchronization process every time you change the batteries.

## Deactivating radio transmission

---

If you want to use or carry your toothbrush in surroundings where radio controlled appliances are not allowed (e.g. hospitals or aircraft), deactivate the radio transmission by simultaneously pressing the on/off and mode button for 3 seconds. The radio transmission light (e) is off while brushing. If you want to reactivate radio transmission, press both buttons again for 3 seconds. The radio transmission light is on while brushing.

Subject to change without notice.

This appliance complies with the requirements covered by the European Directive 1999/5/EC. Receiver Class 3, Duty Cycle Class 3. To find the complete EC declaration of conformity, please visit [www.oralb.com/ce](http://www.oralb.com/ce).

## Environmental Notice

---

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Center or at approved recycling or disposing locations in accordance with local regulations.

### Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects

that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

**UK:**  
FREEPOST RSHJ-KRSA-UHGE  
Braun Product Service  
Olympic Shaver Centre Ltd  
PO BOX 9799  
Calverton  
Nottingham  
NG14 6WN

**ROI:**  
Braun Product Service  
Olympic Shaver Centre Ltd  
Dept AA1758  
PO BOX 6151  
FREEPOST  
Dublin 2

## 30-Days Money Back Guarantee

---

Try Oral-B Pulsonic SmartSeries starting from the day of purchase. If you are not satisfied with this product, please return your handle, Refills, SmartGuide and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

**Great Britain / Ireland:**  
«Oral-B 30-Day Challenge»;  
Dept. PX415  
PO Box 10  
Diss  
Norfolk  
IP98 3SP

For further details, please contact our customer service.

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