
# SWITEL

## BH 700

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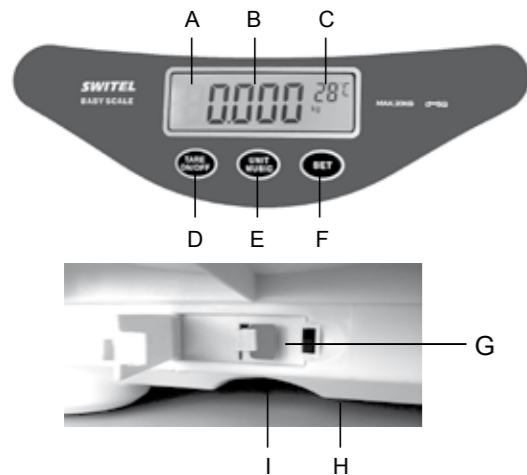#### DEUTSCH

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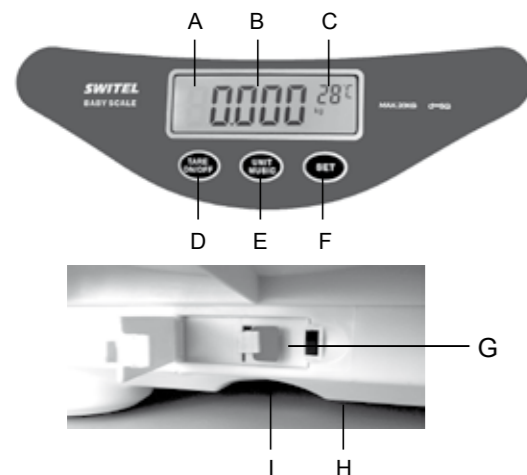
#### ENGLISH

##### Operating Instructions

 Please read and observe the following information and keep the operating manual nearby for future reference!

#### WARNING!

##### Safety information

- ➔ Only use batteries of the same type!
- ➔ Always treat batteries with due care and attention and only use batteries as described.
- ➔ Never use standard and rechargeable batteries together. Never use rechargeable batteries or standard batteries with different capacities or charge statuses together. Never use standard or rechargeable batteries which are damaged.
- Risk of suffocation through small objects, packaging and protective foil!
- ➔ Keep children away from the product and its packaging!
- Risks to health and the environment from batteries!
- ➔ Never open, damage or swallow batteries or allow them to pollute the environment. They could contain toxic and ecologically harmful heavy metals.

#### ATTENTION!

- ➔ Do not place anything with a weight exceeding 30 kg on the scale.
- ➔ Prevent exposure to environmental influences, e.g. smoke, dust, vibration, chemicals, moisture, heat and direct sunlight.
- ➔ Only allow repairs to be completed by qualified personnel.

##### PACKAGE CONTENTS

1 Baby scale, 1 Operating manual

##### DISPLAY ICONS AND OPERATING ELEMENTS

- A Display with backlight
- B Weight (min. 0,1 kg, max. 20 kg, division 5 g, „0-Ld“ = Overload)
- C Temperature
- D TARE / ON / OFF (Turn off = press 3 sec.)
- E UNIT / MUSIC (Press and hold to change from „kg“ to „lb“)
- F SET
- G Measuring tape (on the left side of the unit, max. 150 cm)
- H Battery cover (in the bottom)
- I Reset button for tape measure

##### STARTING UP

- ➔ After inserting the batteries, place the scale on a flat surface.  
**If the scale is not in a stable situation when turned on, the SET-function will not work, the display shows “UnST”**

- ➔ Press **TARE / ON / OFF**. The display shows 0,000 kg
- ➔ Place the baby in the middle of the scale. When the weight value is locked, the music will start to play. Now you can pick up the baby. (After 60 seconds without any operation the scale will automatically switch off.)

##### MEASUREMENT HISTORY

- ➔ Press **SET**  
The Display shows the last weighing operation (P0)
- ➔ Press **SET** for approx. 3 seconds.  
The Display shows the last 10 weighing operations (P1), Press **SET** (P2), Press **SET** (P3), ➔➔➔ Press **SET** (P0)

##### TARE FUNCTION

Place a blanket (max. 10 kg) on the scale if the scale is off.

- ➔ Press **TARE / ON / OFF**. The display shows 0,000 kg

Place a blanket on the scale if the scale is on. The display shows the weight of the blanket.

- ➔ Press **TARE / ON / OFF**. The display shows 0,000 kg

**Take care that the blanket does not touch the ground, as this could disturb the weighing accuracy!**

##### MUSIC PLAYER

###### Change the music

During the weighing process you will hear music playing.

- ➔ Press **UNIT** to change the songs.  
There are 35 songs for your choice.

##### On / Off

**The scale must be off!**

- ➔ Press and hold **SET** and **TARE / ON / OFF**. The display shows „ON or OFF“.
- ➔ Press **UNIT** to change.

##### BATTERIES

If the display shows “Lo” after starting, the batteries need to be replaced.

##### Technical data (All rights reserved)

Feature	Value
BH 700	
W x H x D	650 x 410 x 90 mm
Weight	1.965 g (without battery)
Power	4 x 1,5V AA, not included
Temperature range	-10 ~ +50 °C
Humidity range	to 90 %

##### Service hotline

In the case of technical problems, contact our Service hotline. Switzerland: Tel. 0900 00 1675 (national charges, Swisscom at time of going to print: CHF 2.60/min). In the case of claims under the terms of guarantee, contact your sales outlet.

##### Maintenance

Clean the surface of the housing with a soft, fluff-free cloth. Do not use any cleaning agents or solvents.

##### Guarantee

SWITEL equipment is produced and tested according to the latest production methods. The implementation of carefully chosen materials and highly developed technologies ensure trouble-free functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase. All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair defective parts, replace them or replace the entire device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer. If your device does show signs of a defect within the period of guarantee, please contact the sales outlet where you purchased the SWITEL device, producing the purchase receipt as evidence. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

##### Disposal

 In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to the laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old electronic and electrical devices in a separate waste container. The adjacent symbol indicates that the device must not be disposed of in normal domestic waste!

 You are legally obliged to dispose of power packs and batteries at the point of sale or in the corresponding containers provided at collection points by local public waste authorities. Disposal is free of charge. The symbols indicate that the batteries must not be disposed of in normal domestic waste and that they must be brought to collection points provided by local public waste authorities. Packaging materials must be disposed of according to local regulations.

##### Declaration of Conformity

 This device fulfils the requirements stipulated in the Directive 1999/5/EC on radio equipment and telecommunications terminal equipment and the mutual recognition of their conformity. Conformity with the above mentioned directive is confirmed by the CE mark on the device. To view the complete Declaration of Conformity, please refer to the free download available on our web site [www.switel.com](http://www.switel.com).