

Silk'n™
Beautiful Technology

SatinGlow
H330I

- EN INSTRUCTIONS FOR USE
- DE BEDIENUNGSANLEITUNG
- FR INSTRUCTIONS D'UTILISATION
- NL GEBRUIKSINSTRUCTIES
- ES INSTRUCCIONES DE USO
- IT MANUALE OPERATIVO

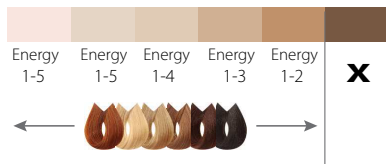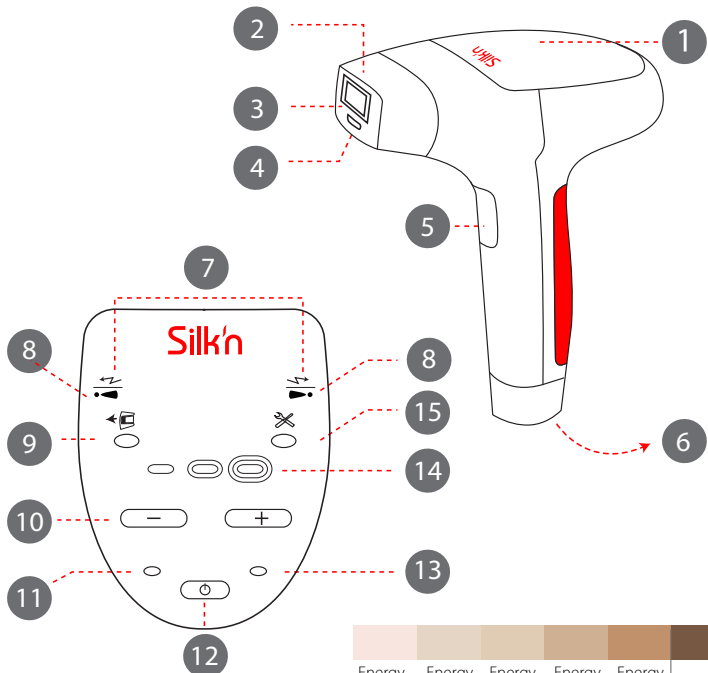Legal Notice

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam POB 533, Yokneam 2069206, ISRAEL.

Table of Contents

1. Before you start.....	3
2. Intended use.....	5
3. Safety with Silk'n SatinGlow	5
4. Contraindications	6
5. Do not use Silk'n SatinGlow if:.....	10
6. Avoid using Silk'n SatinGlow if:.....	11
7. Possible side effects.....	12
8. Long term hair removal with Silk'n	14
9. Get started!.....	16
10. Maintenance.....	23
11. Troubleshooting	25
12. Customer service.....	27
13. Warranty	28
14. Disposal	28
15. Specifications	29
16. Signs and Symbols.....	30

Before using Silk'n SatinGlow for the first time, please read this manual in full. Pay particular attention to sections on device procedures, device operation, and after-use procedures. We recommend you re-familiarize yourself with this manual before each use of Silk'n SatinGlow. Silk'n SatinGlow is a powerful electrical device. As such, it should be used with special attention to safety. Please read all warnings and safety precautions before use, and follow them strictly.

1. Before you start

Silk'n SatinGlow is a light-based device for permanent hair reduction designed for home-use.

1.1 Device description

Your Silk'n SatinGlow consists of a hair removal unit with a *control panel* and a *disposable lamp cartridge*. On the Silk'n SatinGlow hair removal unit handle you can find the *pulse button*. The *disposable lamp cartridge* is located at the tip of the hair removal unit.

1.2 Package contents

Your package contains:

- Silk'n SatinGlow hair removal unit
- Adapter
- User manual

Refer to the diagram at the front of this manual to see the location of each feature.

- | | |
|-----------------------------|---|
| 1 Control panel | 8 Skin tone warning indicator lights |
| 2 Disposable lamp cartridge | 9 Replacement cartridge indicator light |
| 3 Treatment surface | 10 Energy level setting buttons |
| 4 Skin colour sensor | 11 Power on indicator light |
| 5 Pulse button | 12 On/Standby button |
| 6 Adapter cable socket | 13 Standby indicator light |
| 7 Ready indicator light | 14 Energy level indicator lights |
| | 15 System warning indicator light |

1.3 Disposable lamp cartridge

The Silk'n SatinGlow all-purpose *disposable lamp cartridge* can fire a limited and pre-set number of light pulses (regardless of the energy level of these pulses).

For your information a complete series of full-body treatments typically requires 7,000-10,000 light pulses.

The type of *disposable lamp cartridge* included with your system, and the number of light pulses in them, is specified on the package of your Silk'n SatinGlow.

2. Intended use

Silk'n SatinGlow is intended for removal of unwanted hair. Silk'n SatinGlow is also intended to effect long term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs re-growing after a treatment regime.

Silk'n SatinGlow may be used to remove unwanted body hair.

Ideal body areas for Silk'n SatinGlow use include the underarms, bikini line, arms, legs, back shoulders and chest.

3. Safety with Silk'n SatinGlow

With Silk'n SatinGlow safety comes first

Home Pulsed Light™ technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy used in Silk'n SatinGlow reduces its potential to cause harm or complications, and contributes to your overall safety.

Silk'n SatinGlow protects your skin

Light based hair removal is not suitable for naturally darker skin tones or tanned skin. Silk'n SatinGlow comes with a built-in *skin colour sensor* that is designed to measure the complexion of the applied surface and enable application only on suitable skin complexions. This unique safety feature will not let you treat where your skin is too dark or too tanned. Furthermore, the *treatment surface* through which pulses of light are delivered is recessed inside the *lamp cartridge*. This enables Silk'n SatinGlow to protect your skin by avoiding direct contact between the *treatment surface* and the skin.

Silk'n SatinGlow protects your eyes

Silk'n SatinGlow has a built-in safety feature for eye protection. It has been designed so that a light pulse can not be emitted when the *treatment surface* is facing open air. The safety button enables treatment only when the *treatment surface* is in contact with the tissue.

4. Contraindications

Silk'n SatinGlow is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice you may also visit www.silk'n.eu or your local Silk'n domain.

If you are unsure if Silk'n SatinGlow is safe for you to use, please ask your doctor or dermatologist!

Warning	Possible Result	Safety Tip
<p>Do not use Silk'n SatinGlow on naturally dark skin complexion.</p> <p>Not Safe!</p>  <p>See the skin colour chart on the inner cover.</p>	<p>Treating dark skin with Silk'n SatinGlow can result in adverse effects such as redness and skin discomfort.</p>	<p>A unique <i>skin colour sensor</i> is embedded in Silk'n SatinGlow to measure the treated skin complexion at the beginning of each session and occasionally during the session. The <i>skin colour sensor</i> ensures that pulses will only be emitted on suitable skin tones.</p>

Warning	Possible Result	Safety Tip
<p>Do not use Silk'n SatinGlow on tanned skin or after recent sun exposure (including artificial tanning machines). Do not expose treated areas to the sun.</p>	<p>Treating with Silk'n SatinGlow before or after any sun exposure can result in adverse effects such as redness and skin discomfort. Tanned skin particularly following sun exposure, contains large quantities of the pigment Melanin. This applies to all skin types and complexions, including those which don't seem to tan quickly. The presence of large quantities of Melanin exposes the skin to higher risk when using any method of light based hair removal.</p>	<p>Avoid exposure to the sun for 4 weeks before your Silk'n SatinGlow treatment! Sun exposure includes unprotected exposure to direct sunlight of over 15 minutes constantly, or unprotected exposure to diffused sunlight of over 1 hour constantly. To protect recently treated skin when exposed to sunlight, be sure to thoroughly apply sunscreen <i>SPF</i> 30 or higher, for 2 weeks after each hair removal session.</p>
<p>Never use Silk'n SatinGlow around or near the eyes.</p>	<p>Although, Silk'n SatinGlow can be used for the treatment of facial hair, extra caution should be used on the face to avoid the eyes.</p> <p>Possible hazardous optical radiation is emitted from this device. Do not look at the operating lamp. Eye injury may result.</p>	<p>Silk'n SatinGlow has a built-in safety feature for eye protection. It has been designed so that a light pulse can not be emitted when the <i>treatment surface</i> is facing open air. The safety button enables treatment only when the <i>treatment surface</i> is in contact with the tissue.</p>

Warning	Possible Result	Safety Tip
Do not treat on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, birth marks, moles or warts), nipples, genitals or lips.	Treating any area with the Silk'n SatinGlow that is dark in colour or has more pigment can result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypopigmentation).	Cover area with material that will not absorb the light such as a white cloth or white medical tape.
Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).	Light based treatments can cause adverse reactions on already sensitive areas.	Wait for the effected area to heal before using Silk'n SatinGlow.
A treatment area that has a history of herpes outbreaks.	Light based treatments can cause adverse reactions on already sensitive areas.	Consult your physician to received preventative treatment before using Silk'n SatinGlow.
Do not treat with Silk'n SatinGlow if you are pregnant or nursing.	Hormonal changes during pregnancy could increase skin sensitivity and the risk of adverse effects.	
Do not treat on areas where you may want hair to grow.	Results are permanent.	

Warning	Possible Result	Safety Tip
<p>Keep Silk'n SatinGlow away from water! Do not place or store Silk'n SatinGlow where it can fall, be pushed or placed into a tub, sink or any other vessel containing water.</p>	<p>This may cause severe electrocution.</p>	<p>Do not use Silk'n SatinGlow while bathing or if it becomes damp or wet. Unplug immediately if it has fallen into water.</p>
<p>Do not attempt to open or repair your Silk'n SatinGlow device.</p>	<p>Opening Silk'n SatinGlow may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury. It may also damage your device and will void your warranty.</p>	<p>Please contact Silk'n Customer Service if you have a broken or damaged device in need of repair.</p>

5. Do not use Silk'n SatinGlow if:

the <i>disposable lamp cartridge</i> is damaged, or if its <i>treatment surface</i> is cracked, chipped or missing.		Stop use immediately and contact Silk'n Customer Service.
the attachments or accessories are not recommended by Home Skinovations Ltd.		
it has a damaged cord or plug.	Keep the power cord away from heated surfaces.	
you see or smell smoke when it is in use.		
it is not working properly or if it appears to be damaged.		
its <i>skin colour sensor</i> is cracked or coming apart.		
the outer shell is cracked or coming apart.		

6. Avoid using Silk'n SatinGlow if:

you have a history of keloidal scar formation.	Consult with your physician prior to using the Silk'n SatinGlow.
you suffer from epilepsy.	
you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.	
you have a history of skin cancer or areas of potential skin malignancies.	
you have received radiation therapy or chemotherapy treatments within the past 3 months.	
you have any other condition which in your physician's opinion would make it unsafe for you to be treated.	
your skin was treated with a physician-prescribed drug within the past 6 months, please consult with your physician.	

7. Possible side effects

When used according to the instructions, side effects and complications associated with use of Silk'n SatinGlow are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

Adverse Reaction

Minor skin discomfort

Although home pulsed light hair removal is generally very well-tolerated, most users do feel some mild discomfort during use, usually described as being a mild stinging sensation on the treated skin areas. The stinging sensation usually lasts during the time of the application itself or for a few minutes thereafter.

Anything beyond this minor discomfort is and means that either you should not continue to use Silk'n SatinGlow because you are unable to tolerate the hair removal application, or that the energy level setting is too high.

Skin redness

Your skin may become red right after using Silk'n SatinGlow or within 24 hours of using Silk'n SatinGlow. Redness generally clears up within 24 hours. See your doctor if redness does not go away within 2 to 3 days.

Increased skin sensitivity

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Adverse Reaction

Pigment changes

Silk'n SatinGlow targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is risk of temporary hyperpigmentation (increased pigment or brown discolouration) or hypopigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin tones. Usually discolouration or changes to skin pigment are temporary and permanent hyperpigmentation or hypopigmentation rarely occur.

Excessive redness and swelling

In rare cases treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but one should avoid exposure to sun.

Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage treated on the face and neck.

8. Long term hair removal with Silk'n

A typical full hair growth cycle may take 18-24 months. During this time multiple Silk'n SatinGlow sessions may be required in order to achieve permanent hair removal.

The efficiency of hair removal varies from person to person according to body area, hair colour, and how Silk'n SatinGlow is used.

Typical Silk'n SatinGlow hair removal plan during a full hair growth cycle:

- **Treatments 1-4:** Plan two weeks apart
- **Treatments 5-7:** Plan four weeks apart
- **Treatments 8 +:** Treat as needed, until desired results are achieved.

Expected Treatment Outcome for Hair Removal

Fitzpatrick skin type	Hair Color	Light Fluence [J/cm ²]	Anatomic area of the body ²	Average Number of Treatments ¹	Effectiveness of Hair Removal ¹
I-IV	Light brown to brown	2.5 - 4	Lower or upper legs	Session 10-12, approx. four weeks apart.	60%
			Arms		50%
			Bikini line		60%
			Underarm		55%

Fitzpatrick skin type	Hair Color	Light Fluence [J/cm ²]	Anatomic area of the body ²	Average Number of Treatments ¹	Effectiveness of Hair Removal ¹
I-IV	Brown to black	2.5 - 4	Lower or upper legs	Session 8-10, approx. four weeks apart.	70%
			Arms		60%
			Bikini line		70%
			Underarm		65%

1. Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

2. One cannot expect permanent epilation in a single or even double treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

9. Get started!

1. Remove the Silk'n SatinGlow hair removal unit and other components from box.
2. Verify that the *disposable lamp cartridge* is inserted properly at the tip of the device and that it is not damaged.
3. Plug the adapter cord into the Silk'n SatinGlow hair removal unit socket.
4. Plug the adapter into an electrical outlet. The *standby indicator light* will turn on.

Your Silk'n SatinGlow is now ready to start.

-  Treat with Silk'n SatinGlow safely and properly!
 -  Choose areas to treat that have not had recent exposure to the sun.
5. The skin of the treatment areas should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. Do not wax, pluck or tweeze the hair out. Shaving is an important step to get the results you want!
 6. Press the *On/Stdby button*. The *power on indicator light* (orange) will turn on and a fan sound (similar to the sound of a hairdryer) will start.
 7. Approximately 3.5 seconds after pressing the *On/Stdby button*, the *ready indicator light* (green) will turn on. The device is then ready for you to trigger the first pulse at the lowest energy level.

8. If this is your first treatment - Choose the energy level. Silk'n SatinGlow offers 5 levels of energy from 1 (the lowest) to 5 (the highest). To choose the right energy level for your first treatment, test each treatment area separately.
 - » Place the *treatment surface* on the skin and press the *pulse button* to deliver one pulse. If you did not experience abnormal discomfort press the *+ energy level setting button* once and deliver one pulse on a different spot at energy level 2.
 - » Continue testing at increased energy levels until you have reached the highest level that you feel comfortable with. If within one hour you experience no adverse effects, you can start your full treatment at that energy level.
 - » Repeat this test for each body part you intend to treat.
9. Set the energy level. If the energy level should be higher than the lowest, press the *energy level setting buttons* using the "-" or "+" to respectively decrease or increase the energy level, until the desired energy level is set and represented by 3 *energy level indicator lights* (refer to section "9.1 Energy levels").
10. Apply the *treatment surface* to the skin, making sure the skin is spread evenly and smoothly. As soon as the *treatment surface* is in full contact with the skin, the *ready indicator light* will start blinking.
11. Press the *pulse button*. The device will first determine the colour of your skin. If the colour of the skin is light enough for safe application, the device will flash a pulse of light onto your skin, and the *ready indicator light* will turn off.
 - » You will see a bright flash of light and simultaneously hear a subtle pop sound, which is a normal noise for the device. You will feel a mild sensation of warmth and tingling. Silk'n SatinGlow will immediately recharge for the next pulse.
 - » After 3.5 seconds the *ready indicator light* will turn on again.

12. Remove the *treatment surface* from the treated area of skin.
13. Move the *treatment surface* to another area of skin. Use the pressure marks the device just made on your skin to guide you for proper positioning of the next pulse, avoiding both gaps and overlaps between pulses.

Silk'n SatinGlow pulses should be administered in rows, starting at one end of each row and progressing sequentially towards the other end. This technique allows better control of skin coverage, and helps you avoid treating the same area more than once or overlapping skin areas.

Try to avoid overlapping pulses.

Do not treat the same area of skin more than once per hair removal session.

14. Repeat the process until the entire treatment area is covered.

If your skin blisters or burns, stop use immediately!

9.1 Energy levels

The energy level determines the intensity of the Silk'n SatinGlow light pulse delivered to your skin, from the lowest level (-) to the highest level (+). As energy level increases, so does hair removal results as well as the risk of possible side effects and complications.

Silk'n SatinGlow offers 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 3 *energy level indicator lights*:

Level 1 (lowest): The left indicator light is steadily lit.

Level 2: The left indicator light is steadily lit and the central indicator light is blinking.

- **Level 3:** Both left and central indicator lights are steadily lit.
- **Level 4:** The left and central indicator lights are steadily lit and the right indicator light is blinking.
- **Level 5 (highest):** All 3 indicator lights are steadily lit.

Whenever Silk'n SatinGlow is turned on its energy level will automatically be reset to the lowest energy level. Only one *energy level indicator light* will be on.

To set the energy level, press the *energy level setting buttons* using “-” or “+” to respectively decrease or increase the energy level.


9.2 Skin colour sensor

Light-based hair removal on darker skin tones can result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypopigmentation).

To prevent such misuse, a unique *skin colour sensor* in Silk'n SatinGlow measures the treated skin complexion at the beginning of each session and occasionally during the session. If the *skin colour sensor* detects a skin tone that is too dark for Silk'n SatinGlow application, the device will automatically stop emitting pulses.

If you see no light pulse and the *skin tone warning indicator light* blinks, this is an indication that your skin tone, measured by the *skin colour sensor*, is too dark for safe application. Try using the device on a different body part, or contact Customer Service.

9.3 Pulsing and gliding

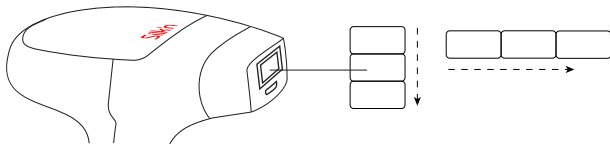Pulsing is the best technique for working at higher energy levels.

- Press the *treatment surface* against your skin. Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end.
- The *treatment surface* creates temporary pressure marks on the treated area so you can see areas you have treated.
- After each pulse, move the *treatment surface* to the next spot. Keep doing this pulsing action for the area you are treating.
- Avoid gaps or overlaps as you use the device.

Gliding is the best technique for working at lower energy levels.

- Press the *treatment surface* firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the *treatment surface* is always touching your skin.
- At lower energy levels, the device takes only 1-3.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.
- Glide at a speed so that you avoid gaps or overlaps.

- ! Do not treat the same area of skin more than once per hair removal session.
- ! Stop use immediately if you have any discomfort.

9.4 What to expect with Silk'n SatinGlow

For many people, using Silk'n SatinGlow may be their first experience with a light-based device designed for home use. Silk'n SatinGlow is simple to use, and hair removal sessions go by quickly. During a Silk'n SatinGlow session it is normal to experience and feel:

- **A fan noise:** The cooling fan in Silk'n SatinGlow makes a noise similar to a hairdryer. This is normal.
- **A pop sound with each pulse:** When a pulse of HPL™ light is activated, it is normal to hear a subtle pop sound simultaneously with the flash of light.
- **Moderate pressure of the treatment surface:** This is necessary and helpful for placement of adjacent pulses of light, and is part of the unique safety feature of Silk'n SatinGlow.
- **A sensation of warmth and tingling:** During each pulse of light it is normal to feel a mild sensation of warmth and tingling from the light energy.
- **Some mild red or pink color:** During and just after your Silk'n SatinGlow session it is not uncommon to see some very mild, pink-like colour of the skin. This is usually most noticeable around the hairs themselves.
- ! However if you see full redness of the skin, blistering or burns, stop using Silk'n SatinGlow immediately.

9.5 After treating with Silk'n SatinGlow

- When you have completed your treatment session, turn off Silk'n SatinGlow by pressing the *On/Stop* button. All the lit indicator lights will turn off and the *Stop* indicator light will turn on. Be sure to remember the last energy level setting you used, as it will not be restored when turning Silk'n SatinGlow on again.
- Unplug the adapter cord from the electrical outlet.
- After each hair removal session it is recommended that you clean your Silk'n SatinGlow device, especially the *treatment surface* (refer to section "10.1 Cleaning").
- After cleaning, it is recommended to store your device in its original box, and keep it away from water.

Skin care following a hair removal session:

- ☀ Do not expose treated areas of skin to the sun. Be sure to carefully protect the treated skin with sunscreen, throughout the hair removal period and for at least 2 weeks following the last Silk'n SatinGlow session.

10. Maintenance

10.1 Cleaning

After each hair removal session, it is recommended to clean your Silk'n SatinGlow device, and especially the *treatment surface*.

- Unplug Silk'n SatinGlow before cleaning.
- Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe Silk'n SatinGlow surface, and especially the *treatment surface*.
- ! Never immerse Silk'n SatinGlow or any of its parts in water!

10.2 Replacing the lamp cartridge

The Silk'n SatinGlow *disposable lamp cartridge* can fire a pre-set and limited number of light pulses. Pulse intensity is determined only according to the energy level setting of the device. There is no decrease of energy during the usable lifetime of the *lamp cartridge*.

- When only 100 light pulses remain in the *disposable lamp cartridge*, the *replacement cartridge indicator light* on the *control panel* will start blinking, indicating that the *disposable lamp cartridge* should soon be replaced.
- As the number of light pulses in the *disposable lamp cartridge* further decreases, the *replacement cartridge indicator light* will blink more quickly.
- Once all the pulses in a *disposable lamp cartridge* have been used, the *replacement cartridge indicator light* will remain steadily on, indicating that further pulses cannot be emitted.

- In order to continue the hair removal session, it will now be necessary to replace the *disposable lamp cartridge*.

 A *disposable lamp cartridge* should also be replaced if large spots appear inside it or if its *treatment surface* is broken. For replacement parts, contact Customer Service.

To replace a *disposable lamp cartridge*:

1. Press the *On/Stdby button* to turn Silk'n SatinGlow to standby mode.
2. Unplug the power adapter from the electrical outlet.
3. Grasp the used *disposable lamp cartridge* on both sides, pull it out of the socket and discard with normal trash.
4. Unwrap a new *disposable lamp cartridge*.
5. Push the new *disposable lamp cartridge* gently into place.

 After replacing the *disposable lamp cartridge* with a new one, always set the energy to be one level lower than your normal setting. Raise the energy level back, only if you do not experience abnormal discomfort during or after the treatment.

11. Troubleshooting

11.1 “My Silk’n SatinGlow does not start.”

- Make sure the power adapter is properly connected to the Silk’n SatinGlow device.
- Make sure the power adapter is plugged into an electrical outlet on the wall.

11.2 “A light pulse is not emitted when I press the Pulse Button.”

Make sure that the *ready indicator light* is blinking to ensure you have good contact with the skin and that the *treatment surface* is evenly and firmly pressed to the skin. For your safety, the *pulse button* will activate a pulse only if the *treatment surface* is firmly pressed against the skin.

- Check the *skin tone warning indicator light*. If it blinks, this is an indication that your skin tone, measured by the *skin colour sensor*, is too dark for safe application. Try using the device on a different body part or contact Silk’n support.
- Check the *replacement cartridge indicator light*. If it is on, disconnect Silk’n SatinGlow and replace the *disposable lamp cartridge*.
- Make sure that the *ready indicator light* is on.
 - » If within 10 seconds the *ready indicator light* remains off, turn Silk’n SatinGlow off and back on by pressing the *On/Stdby button* twice.
 - » If the problem persists, contact your local Silk’n Customer Service Center.

- Check the *system warning indicator light*.
 - » If the light is on, turn Silk'n SatinGlow off and back on by pressing the *On/Stdby button* twice.
 - » If the problem persists, contact your local Silk'n Customer Service Centre.
- ! Do not attempt to open or repair your device. Only authorized Silk'n repair centers are permitted to perform repairs.

Opening Silk'n SatinGlow may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury.

Trying to open Silk'n SatinGlow may also damage the device and will void your warranty. Please contact Silk'n Customer Service if you have a broken or damaged device in need or repair.

12. Customer service

See our website www.silkn.eu for more information. This manual is also available as a PDF download from the website. Contact Silk'n Customer Service to report unexpected operations or events, or any other problem with your device.

Our contact details are:

Country	Service number	Email
Belgium	0900-25006	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservice@inno-essentials.de
Italy		servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicenl@silkn.eu
Spain	807088200	servicioalcliente@silkn.eu
UK	0906-2130009	customerareuk@silkn.eu
US and Canada	1-877-367-4556	contact@silkn.com
Other countries	+31 (0)180-330 550	info@silkn.com

13. Warranty

Refer to your Silk'n warranty card.

14. Disposal

Help protect the environment. Do not dispose with household waste and be sure to follow your country's laws for the disposal of electrical items.

15. Specifications

Spot size	1.3cm x 3cm [3.9cm ²]	
Speed	1 pulse every 3.5 second: 1.1 cm ² /sec	
Technology	Home Pulsed Light™	
Max Energy Level	Max 4J/cm ²	
Wavelength	475-1200nm	
Charging Time/Power Source	Continuous operation	
Electrical Requirements	100-240 VAC, 0.9A	
Operation and Safety	Skin colour sensor ensures use on suitable skin types	
Package Size	22cm x 22cm x 13.5cm	
System weight	0.35kg	
Temperature	Operating	10 to 40°C
	Storage	-40 to +70°C
	Operating	30 to 75% rH
Relative Humidity	Storage	10 to 90% rH
Atmospheric pressure	Operating	700 to 1060 hPa
	Storage	500 to 1060 hPa

16. Signs and Symbols

	Products marked with this symbol fulfil the requirements of the EU directives.
	Class II. The adapter is protected against electric shock by means of double insulation.
IP30	The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.
	Do not use in wet environments.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.

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Copyright 2017 © Home Skinovations Ltd.